

FCR100

Intelligent Stepper Motor Rotation Stages



Warranty

Newport Corporation warrants this product to be free from defects in material and workmanship for a period of 1 year from the date of shipment. If found to be defective during the warranty period, the product will either be repaired or replaced at Newport's discretion.

To exercise this warranty, write or call your local Newport representative, or contact Newport headquarters in Irvine, California. You will be given prompt assistance and return instructions. Send the instrument, transportation prepaid, to the indicated service facility. Repairs will be made and the instrument returned, transportation prepaid. Repaired products are warranted for the balance of the original warranty period, or at least 90 days.

Limitation of Warranty

This warranty does not apply to defects resulting from modification or misuse of any product or part.

CAUTION

Warranty does not apply to damages resulting from:

- **Incorrect usage:**
 - Driven load greater than maximum specified load.
 - Stage speed higher than specified.
 - Improper grounding.
 - Connectors must be properly secured.
 - When the load on the stage represents an electrical risk, it must be connected to ground.
 - Excessive or improper cantilever loads.
 - **Modification of the stage or any part thereof.**
-



CAUTION

Please return equipment in the original (or equivalent) packing.

Newport will not be responsible for damage incurred from inadequate packaging if the original packaging is not used.

This warranty is in lieu of all other warranties, expressed or implied, including any implied warranty of merchantability or fitness for a particular use. Newport Corporation shall not be liable for any indirect, special, or consequential damages.

©2018 by Newport Corporation, Irvine, CA. All rights reserved.

Original instructions.

No part of this document may be reproduced or copied without the prior written approval of Newport Corporation. This document is provided for information only, and product specifications are subject to change without notice. Any change will be reflected in future publishings.

Table of Contents

Warranty	ii
EC Declaration of Conformity	v
Definitions and Symbols	vi
Warnings	vii
Cautions	viii
<hr/>	
1.0 Introduction	1
<hr/>	
2.0 System Overview.....	2
2.1 General Description	2
2.1.1 Package Contents.....	2
2.1.2 Controller Specifications	3
2.1.3 Design Details	3
2.1.4 Dimensions.....	4
2.1.5 Weight	4
2.1.6 Assemblies.....	5
2.1.7 FC-PS40 Power Supply	6
2.1.8 USB-RS422-1.8 USB Adapter	6
2.1.9 FC-CB1 Daisy Chain Communication Cable.....	6
2.1.10 EQ120 Bracket	7
2.2 Characteristics.....	8
2.2.1 Definitions	8
2.2.2 Mechanical Specifications.....	9
2.2.3 Load Characteristics and Stiffness.....	9
2.3 System Environmental Specifications	9
2.4 Connector Identification.....	10
2.5 Communication Settings	11
<hr/>	
3.0 Getting Started.....	12
3.1 Mounting and Initial Setup.....	12
3.2 Connection	12
3.3 Initialization and Homing.....	14
3.4 Daisy-Chaining	15
3.5 Multi Stage Configuration.....	16
<hr/>	
4.0 Programming	17
4.1 State Diagram.....	17
4.2 Command Syntax.....	19

4.3	Command Execution Time.....	19
4.4	Command Set.....	20
	AC — Set/Get acceleration	22
	BA — Set/Get backlash compensation.....	23
	BH — Set/Get hysteresis compensation.....	24
	FR — Set/Get stepper motor configuration.....	25
	HT — Set/Get HOME search type	26
	ID — Set/Get stage identifier	27
	JR — Set/Get jerk time.....	28
	MM — Enter/Leave DISABLE state.....	29
	OH — Set/Get HOME search velocity.....	30
	OR — Execute HOME search.....	31
	OT — Set/Get HOME search time-out.....	32
	PA — Move absolute	33
	PR — Move relative	34
	PT — Get motion time for a relative move	35
	PW — Enter/Leave CONFIGURATION state.....	36
	RS — Reset controller	37
	RS## — Reset controller's address	38
	SA — Set/Get controller's RS422 address.....	39
	SE — Configure/Execute simultaneous started move	40
	SL — Set/Get negative software limit.....	41
	SR — Set/Get positive software limit.....	42
	ST — Stop motion.....	43
	TB — Get command error string.....	44
	TE — Get last command error.....	45
	TH — Get set-point position	46
	TP — Get current position.....	47
	TS — Get positioner error and controller state.....	48
	VA — Set/Get velocity.....	50
	VE — Get controller revision information.....	51
	ZT — Get all configuration parameters.....	52
<hr/>		
5.0	Maintenance	53
5.1	Maintenance	53
5.2	Repair	53
5.3	Troubleshooting	54
<hr/>		
	Service Form	55

EC Declaration of Conformity

FCL Series


Newport®

EU Declaration of Conformity

following Annex II-1A
of Directive 2006/42/EC on machinery

The manufacturer:

MICRO-CONTROLE Spectra-Physics,
9, rue du bois sauvage
F-91055 Evry FRANCE

Hereby declares that the machinery:

- Description: " FCL "
- Function: Low Cost Linear Stage
- Models: FCL50 / FCL 100 / FCL 200

– the technical file of which was compiled by:
Mr Hervé LE COINTE , Quality Director,
MICRO-CONTROLE Spectra-Physics, Zone Industrielle - B.P.29
F-45340 Beaune La Rolande France

– complies with all the relevant provisions of the Directive 2006/42/EC on machinery.
– complies with all the relevant provisions of the Directive 2014/30/EU relating to electro-magnetic compatibility.
– complies with all the relevant provisions of the Directive 2011/65/EU relating to RoHS2.

– was designed and built in accordance with the following harmonised standards:

- NF EN 61326-1:2013 « Electrical equipment for measurement, control and laboratory use – EMC requirements – Part 1: General requirements »
- NF EN 55011:2010/A1:2011 Class A

– was designed and built in accordance with the following other standards:

- NF EN 61000-4-2
- NF EN 61000-4-3
- NF EN 61000-4-4
- NF EN 61000-4-5
- NF EN 61000-4-6

ORIGINAL DECLARATION

Done in Beaune La Rolande on 16 May 2017
Hervé LE COINTE
Quality Director



DC1-EN rev:A

Definitions and Symbols

The following terms and symbols are used in this documentation and also appear on the product where safety-related issues occur.

General Warning or Caution



The exclamation symbol may appear in warning and caution tables in this document. This symbol designates an area where personal injury or damage to the equipment is possible.

The following are definitions of the Warnings, Cautions and Notes that may be used in this manual to call attention to important information regarding personal safety, safety and preservation of the equipment, or important tips.



WARNING

Warning indicates a potentially dangerous situation which can result in bodily harm or death.



CAUTION

Caution indicates a potentially hazardous situation which can result in damage to product or equipment.

NOTE

Note indicates additional information that must be considered by the user or operator.

European Union CE Mark



The presence of the CE Mark on Newport Corporation equipment means that it has been designed, tested and certified as complying with all applicable European Union (CE) regulations and recommendations.

Warnings and Cautions



ATTENTION

This stage is a Class A device. In a residential environment, this device can cause electromagnetic interference. In this case, suitable measures must be taken by the user.

Warnings



WARNING

The rotation of objects of all types carries potential risks for operators. Ensure the protection of operators by prohibiting access to the dangerous area and by informing the personnel of the potential risks involved.

WARNING

Do not use this stage when its motor is emitting smoke or is unusually hot to the touch or is emitting any unusual odor or noise or is in any other abnormal state.

Stop using the stage immediately, switch off the motor power by disconnecting the power supply.

After checking that smoke is no longer being emitted contact your Newport service facility and request repairs. Never attempt to repair the stage yourself as this can be dangerous.

WARNING

Make sure that this stage is not exposed to moisture and that liquid does not get into the stage.

Nevertheless, if any liquid has entered the stage, switch off the motor power by disconnecting the power supply.

Contact your Newport service facility and request repairs.



WARNING

Do not insert or drop objects into this stage, this may cause an electric shock, or lock the drive.

Do not use this stage if any foreign objects have entered the stage. Switch off the motor power by disconnecting the power supply.

Contact your Newport service facility for repairs.

WARNING

Do not place this stage in unstable locations such as on a wobbly table or sloping surface, where it may fall or tip over and cause injury.

If this stage has been dropped or the case has been damaged, switch off the motor power by disconnecting the power supply.

Contact your Newport service facility and request repairs.

WARNING

Do not attempt to modify this stage; this may cause an electric shock or downgrade its performance.

WARNING

Do not exceed the usable depth indicated for the mounting holes (see section “Dimensions”). Longer screws can damage the mechanics or cause a short-circuit.

Cautions

CAUTION

Do not place this stage in a hostile environment such as X-Rays, hard UV,... or in any vacuum environment.

CAUTION

Do not place this stage in a location affected by dust, oil fumes, steam or high humidity. This may cause an electric shock.

CAUTION

Do not leave this stage in places subject to extremely high temperatures or low temperatures. This may cause an electric shock.

- Operating temperature: +10 to +35 °C.
- Storage temperature: -10 to +40 °C (in its original packaging).



CAUTION

Do not move this stage if its motor power is on.

Make sure that cables are disconnected before moving the stage. Failure to do so may damage the cable and cause an electrical shock.

CAUTION


Be careful that the stage is not bumped when it is being carried. This may cause it to malfunction.

CAUTION

When handling this stage, always unplug the equipment from the power source for safety.

CAUTION

Contact your Newport service facility to request cleaning and specification control every year.

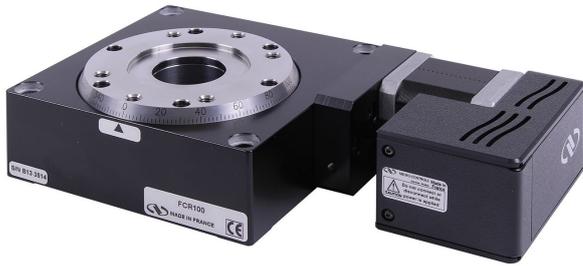

Intelligent Stepper Motor Rotation Stages FCR100

1.0 Introduction

This manual provides operating instructions for the FCR100 intelligent stepper motor rotation stage that you have purchased.

The FCR100 is a precision rotation stage with an integrated stepper motor/controller. The stage can be easily controlled from a PC, using a downloadable graphical user interface (GUI). Communication with the FCR100 is achieved via an RS-422 serial link. A USB to RS422 adaptor can be used (requires Windows™ XP or Windows™ 7 operating system). While the GUI enables basic motion, advanced application programs can be developed via an ASCII command interface and a set of two letter mnemonic commands.



FCR100 rotation stage.

2.0 System Overview

2.1 General Description


The FCR100 is a compact, value-engineered rotation stage with integrated stepper motor and controller. By integrating the electronics into the stepper motor housing, the Newport iPP intelligent stepper motor simplifies the setup greatly, eliminating the need for separate controller, driver box or motor cables. Building on the robust construction of crossed roller bearing and worm screw, the FCR100 provides high motion sensitivity and position accuracy. With the stepper motor directly attached to the worm screw, position is attained by the number of commanded steps or micro-steps without encoder feedback. Up to 4 FC series units can be daisy-chained and controlled by a single GUI, via USB (USB-RS422-1.8) and direct RS-422 (daisy chaining cable FC-CB1).

2.1.1 Package Contents

- FCR100 Intelligent Stepper Motor Rotation Stage
- Stage test report

The following parts are to be ordered separately

- FC-PS40 FC Stage Power Supply (Cable Length: 1.75 m)
- USB-RS422-1.8 USB Adapter (Cable Length: 1.8 m)
- FC-CB1 Daisy Chain Communication Cable (Cable Length: 1 m)

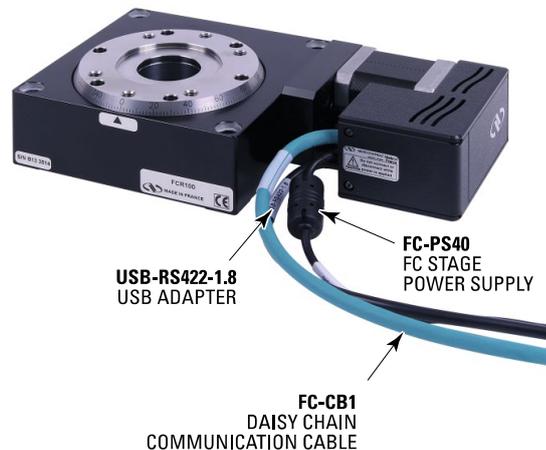


2.1.2 Controller Specifications

General description	Single Axis Intelligent Stepper Motor
Control capability	Stepper motors, open loop
Output power	24 VDC peak, 2.5 A peak PWM amplifier
Control loop	<ul style="list-style-type: none"> – PI control of motor phase current with programmable idle state – 2 kHz servo rate – Backlash compensation
Motion	Point-to-point motion with S-gamma profile and jerk time control.
Computer interface	USB +5 V (USB): < 0.5 A , Windows Compatible
Programming	<ul style="list-style-type: none"> – 40+ intuitive, 2-letter ASCII commands – Command set includes software limits, user units, synchronized motion start, stop all...
Communication rate	115 200 baud (serial link over USB)
Internal safety features	End of range checks, power limit checks, watchdog timer
Consumption	+24 V (FC-PS40): < 1.67 A


2.1.3 Design Details

Base Material	Aluminum
Motor	Integrated iPP 2-phase step motor and controller/driver
Drive Mechanism	Ground worm gear
Worm gear ratio	1:200
Bearing	Crossed roller bearings
Limit Switches	None
Origin	Optical
Communication	RS422, USB
Manual Adjustment	Via 2 mm Allen wrench imprint at the end of the worm screw
Cable	Cables and power supply sold separately

2.1.4 Dimensions

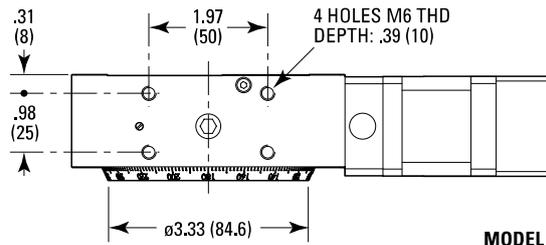

MODEL SHOWN: FCR100
DIMENSIONS IN INCHES (AND MILLIMETERS)

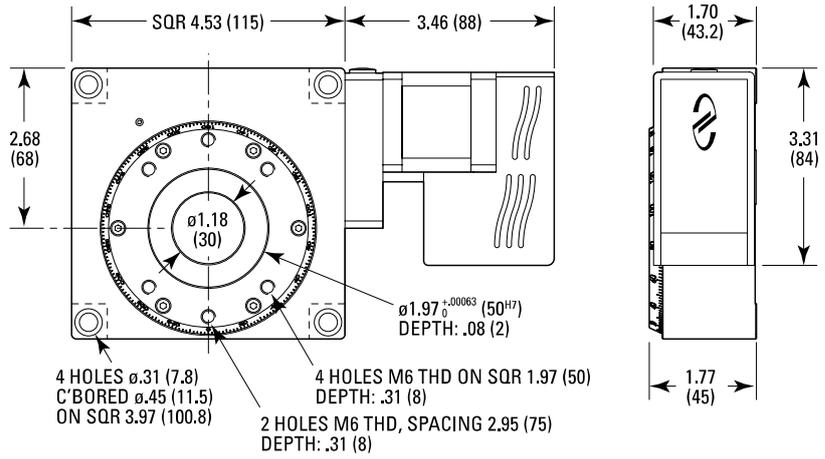

2.1.5 Weight

	Weight [lb (kg)]
FCR100	5.0 (2.25)

2.1.6 Assemblies

Below are a few assembly examples with FC linear and rotation stages. The FC series mounting interfaces use common hole patterns, eliminating the need for adapter plates. Call for compatibility with other Newport stages.



FCL200 and FCR100 rotation stage.



2 FCL stages in XY configuration and one FCR100 rotation stage.



2 FCL stages in XZ configuration with an EQ120 bracket.

2.1.7 FC-PS40 Power Supply

NOTE

Each FC series stage requires a FC-PS40 power supply.



Box Dim. (L x D x H)	125 x 50 x 32 mm
----------------------	------------------

Specifications

AC Input	100–240 VAC, 50–60 Hz, 1 A
DC Output	24 V, 1.67 A, 40 W max.
Cable Length	1.75 m between the power supply box and the iPP driver board connector

2.1.8 USB-RS422-1.8 USB Adapter




Cable Length	1.8 m
--------------	-------

2.1.9 FC-CB1 Daisy Chain Communication Cable



Cable Length	1 m
--------------	-----

2.1.10 EQ120 Bracket

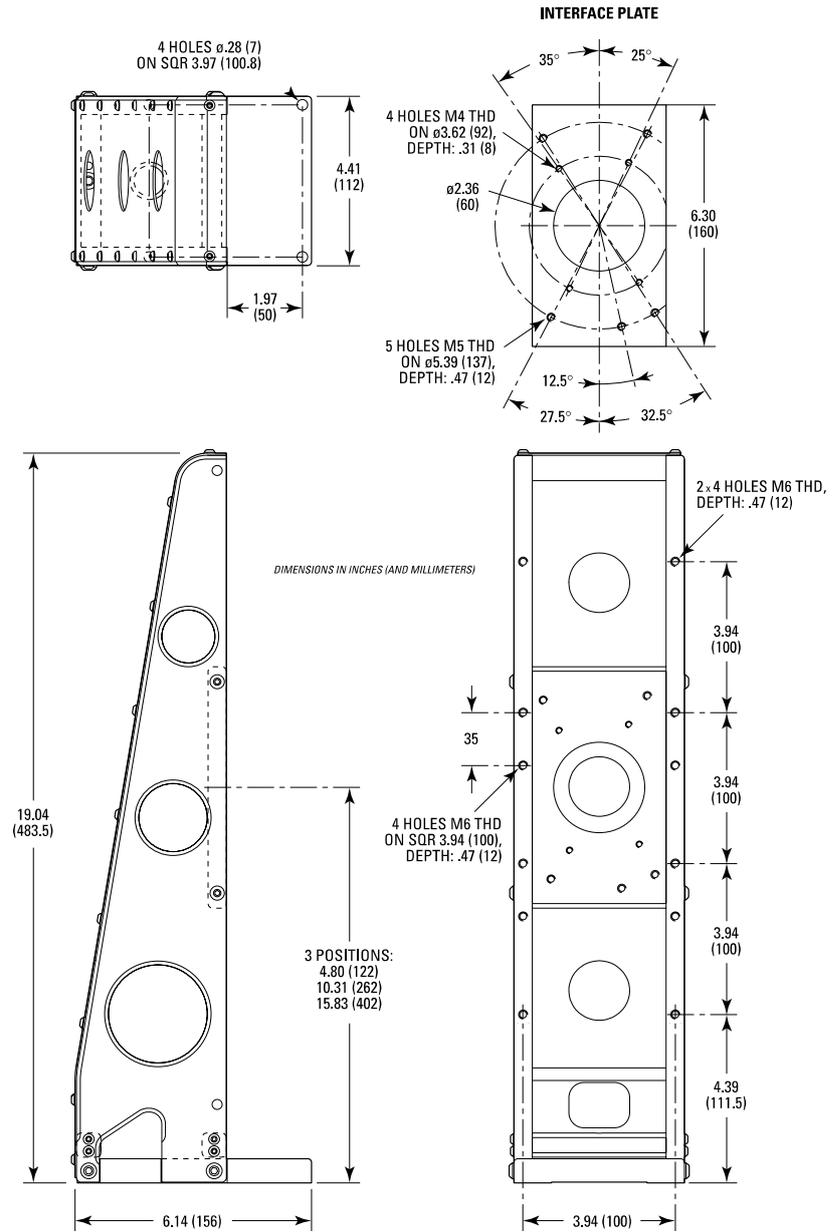

2.2 Characteristics

2.2.1 Definitions

Specifications of our products are established in reference to ISO 230 standard part II “Determination of accuracy and repeatability of positioning numerically controlled axes”.

This standard gives the definition of position uncertainty which depends on the 3 following parameters:

(Absolute) Accuracy

Difference between ideal position and real position.

On-Axis Accuracy

Difference between ideal position and real position after the compensation of linear errors.

Linear errors include: cosine errors, inaccuracy of screw or linear scale pitch, angular deviation at the measuring point (Abbe error) and thermal expansion effects. All Newport motion electronics can compensate for linear errors.

The relation between absolute accuracy and on-axis accuracy is as follows:

$$\text{Absolute Accuracy} = \text{On-Axis Accuracy} + \text{Correction Factor} \times \text{Travel}$$

Repeatability

Ability of a system to achieve a commanded position over many attempts.

Reversal Value (Hysteresis)

Difference between actual position values obtained for a given target position when approached from opposite directions.

Minimum Incremental Motion (MIM or Sensitivity)

The smallest increment of motion a device is capable of delivering consistently and reliably.

Resolution

The smallest increment that a motion device can theoretically move and/or detect. Resolution is not achievable, whereas MIM, is the real output of a motion system.

Eccentricity

Displacement of the geometric center of a rotation stage from the rotation axis in the plane defined by bearings.

Wobble

Tilt of rotation axis during rotation of a stage, measured on a reference surface.

The testing of on-axis accuracy, repeatability, and reversal error are made systematically with test equipment in a controlled room (20 ±1 °C).

The test consists of 4 cycles in each direction, with 21 data points over the travel resulting in a total 164 data points.

Guaranteed Specifications

Guaranteed maximum performance values are verified per Newport's A167 metrology test procedure. For more information, please consult the metrology tutorial section in the Newport catalog or at www.newport.com.

2.2.2 Mechanical Specifications



Travel Range (°)	360 continuous
Minimum Incremental Motion (°)	0.00025
Uni-directional Repeatability (°)	0.002 or ± 0.001
Bi-directional Repeatability, Guaranteed ⁽¹⁾ (°)	0.012 or ± 0.006
Absolute Accuracy, Guaranteed ⁽¹⁾ (°)	0.04 or ± 0.02
Maximum Speed (°/s)	20
Wobble, Guaranteed ⁽¹⁾ (μ rad)	50 or ± 25
Eccentricity, Guaranteed ⁽¹⁾ (μ m)	6 or ± 3
MTBF (h)	20,000


1) Shown are peak to peak, guaranteed specifications or \pm half the value as sometimes shown. The typical specifications are about 2X better than the guaranteed values.



WARNING

To achieve the guaranteed specifications stated, stages must be fixed to a plane surface with a flatness of 5 μ m or better.

2.2.3 Load Characteristics and Stiffness

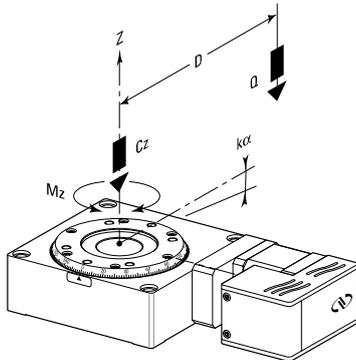


Cz, Normal center load capacity on bearings	300 N
$k\alpha$, Transversal compliance	12 μ rad/Nm
Mz, Maximum torque	0.5 Nm
Jz, Maximum inertia	0.1 kg.m ²
Q, Off-center load	$Q \leq Cz/(1 + D/30)$
D, Cantilever distance in mm	

2.3 System Environmental Specifications

Operating temperature	5 °C to 40 °C
Operating humidity	20% to 85% relative humidity, non-condensing
Location	Indoor use only

2.4 Connector Identification

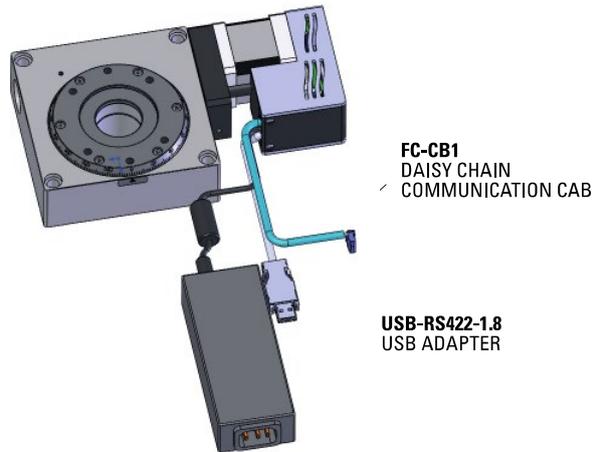


Each FC series stage includes two RS422 communication connectors and a 24 V DC input connector. RS422 connector is connected to either PC (with USB-RS422-1.8 cable) or to RS422 connector of another FC series stage (by FC-CB1). Two RS422 connectors are interchangeable: either of the two connectors can be used for PC connection or the input/output for daisy-chaining connection.



RS422	(2) 6-Pin RS422 communication connectors (connect to USB-RS422-1.8 or to FC-CB1 cable)
24 V in	3-Pin JST PAP-03V-S (connect to FC-PS40)
Cable retainer	2 x M3 threaded hole to attach cable retainer

PINOUT	
1	N.C
2	GND
3	+Tx
4	-Tx
5	-Rx
6	+Rx

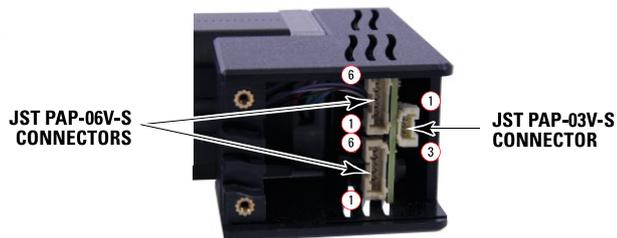

PINOUT	
1	N.C
2	24 V
3	GND

2.5 Communication Settings

Communication parameters are preset in the FC series controller and do not require any configuration:


Bits per second	115 200
Data bits	8
Parity	None
Stop bits	1
Flow control	None
End of line terminator	C _R L _F

- Communication standard: RS-422 4-wire full duplex without handshaking.
- 120 Ω termination resistor on receiver channel.

3.0 Getting Started

3.1 Mounting and Initial Setup

To achieve the guaranteed specifications, stages must be fixed on a plane surface with a flatness of 5 μm .

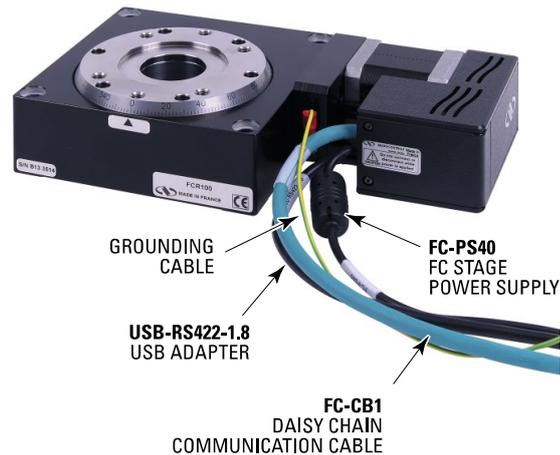

WARNING

The FC series, via threaded hole (M4), must be grounded to avoid electrical disturbances generated by ground loops.

3.2 Connection



WARNING

Do not connect or disconnect cables to FC series stage while power is applied.

Follow the steps described below for the cable connections:


Remove the cover from integrated motor/controller housing by loosening the two screws.



Refer to the Chapter 2.1 and Chapter 2.4 for the description of cables and connectors.

Make proper connections to RS422 connector and 24 VDC power connector, using the USB-RS422-1.8 USB adapter and the FC-PS40 cable respectively.

The second RS422 connector is used only for daisy-chaining.

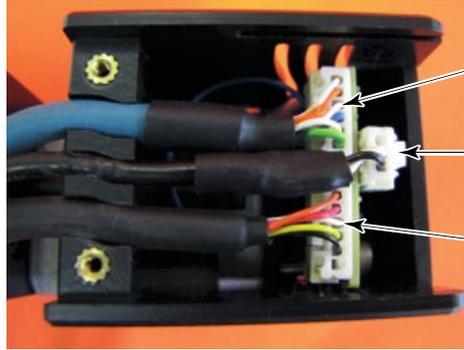


Connect the USB-RS422-1.8 USB to a PC first. Then, connect the FC-PS40 power supply to an electrical outlet.



WARNING

The Power supply cable must remain fully disconnected from electrical outlet when you plug the USB-RS422-1.8 USB adaptor cable to the PC.


When completed, assemble the cover back into integrated motor/controller housing.

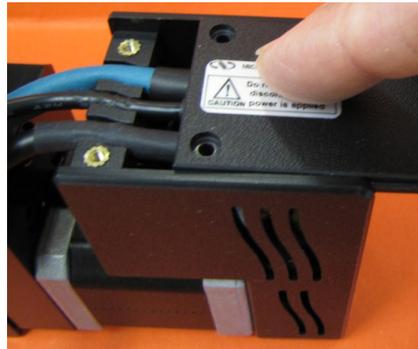

To start using the FCR100, refer to the FC series controller GUI manual.

3.3 Initialization and Homing

In the Initialization and Configuration area, the first button “Initialize” executes homing and changes the controller state from “Not Referenced” to “Ready”. Once it is in the “Ready” state, the “Disable” or “Enable” buttons change the controller state to “Disabled” or “Ready”.

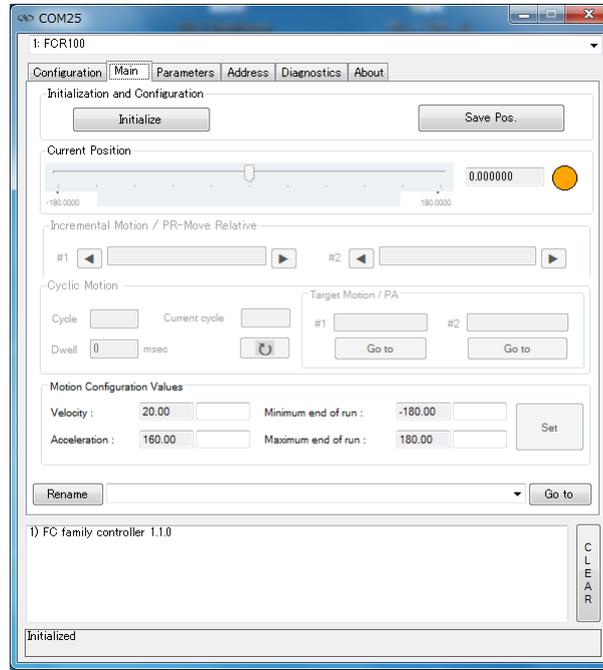

WARNING

The FCR home routine does not behave the same way as other Newport rotation stages. Please read this section.

For the different controller states, refer to section 1.3 of the FC Controller GUI Manual. For further details of the Main tab of the GUI, refer to section 3.2 of the FC Controller GUI Manual.

Due to the unique design of the FCR100, the home process is not the same based on the last position before the homing routine.

- When the last position prior to homing is between -23° and $+180^\circ$, the stage will move directly towards the origin. For example, at -22° , the stage will home in the positive direction (CW) straight to 0° . If the stage is at 125° , the stage will home in the negative direction (CCW) straight to 0° . In both cases, the stage will not cross the negative software limit.
- However, when the last position prior to homing is anywhere from -180° to -23° , homing will move in the negative direction (CCW) and it will cross the negative software limit, until it reaches the origin at 0° . Note: the actual range might be slightly different from $[-180; -23]$ due to the mounting tolerance of origin switch.



WARNING

In cases where interference or collision is an issue beyond the -23° software limit, make sure that the last position of the FCR100 stage before homing is between -23° and $+180^\circ$.


3.4 Daisy-Chaining

Up to 4 FC series controllers can be networked through the internal RS422 communications link. Before daisy-chaining, the controller address of each stage must be set separately via the USB and the GUI. The FC controller that will be directly connected to the PC must have its controller address set to 1 and all subsequent stages must have a different controller address set between 2 and 4. Refer to the GUI manual or command interface manual for instructions to change the controller address.



Once the controller address of each stage is set, unplug all the power cables and disconnect the USB cable from all the stages except for the FC controller that has its address set to 1. Use the daisy chain cables to connect each controller. For stages #2 thru 4, it does not matter whether you use the upper or lower RS422 connector to daisy chain.

- Only the FC controller with address 1 should be connected to the PC
- All other FC controllers must be daisy-chained.

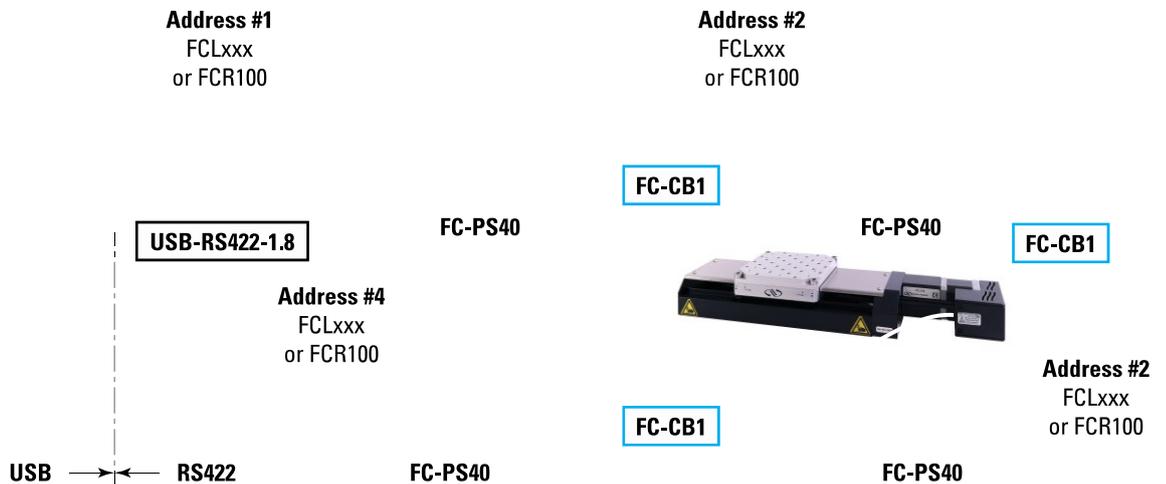

ter


- Each stage from the chain must have a different address from 1 to 4 (default is 1).

NOTE

Another multi-axis option is to run multiple GUI's, each with an FCR100 or daisy-chained FCR100's.

3.5 Multi Stage Configuration

Another way to control multiple stages is opening multiple GUI's for every FC stage. Each of the FC stages has to be discovered independently.

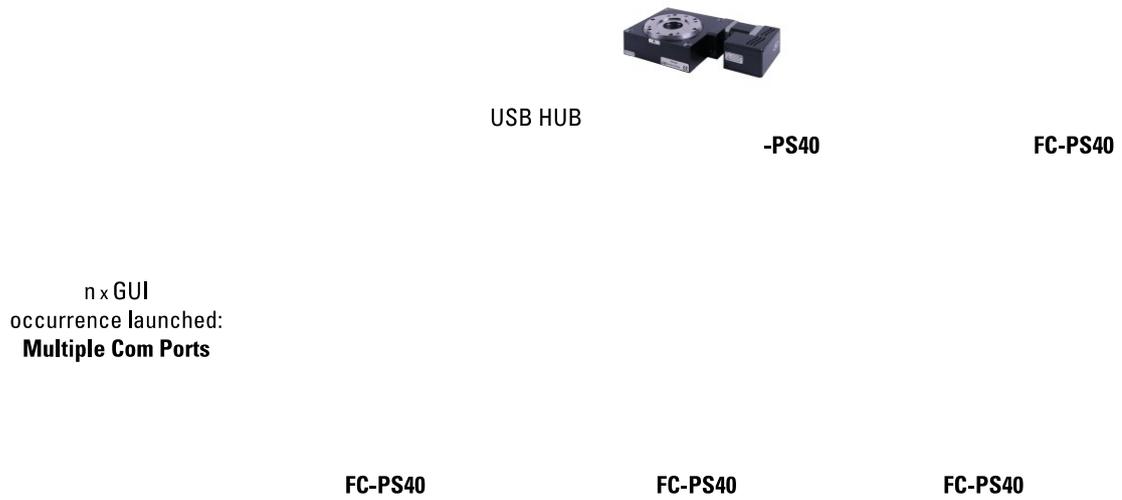


- Address #1 (default) for all stages

4.0 Programming

4.1 State Diagram

For a safe and consistent operation, the FC series stage uses 6 different operational states: Not referenced, Configuration, Homing, Ready, Disable and Moving. In each state, only specific commands are accepted by the FC series stage. Therefore, it is important to understand the state diagram below and to know which commands and actions cause transitions between the different states. Also see section 4.4 for command/state information:

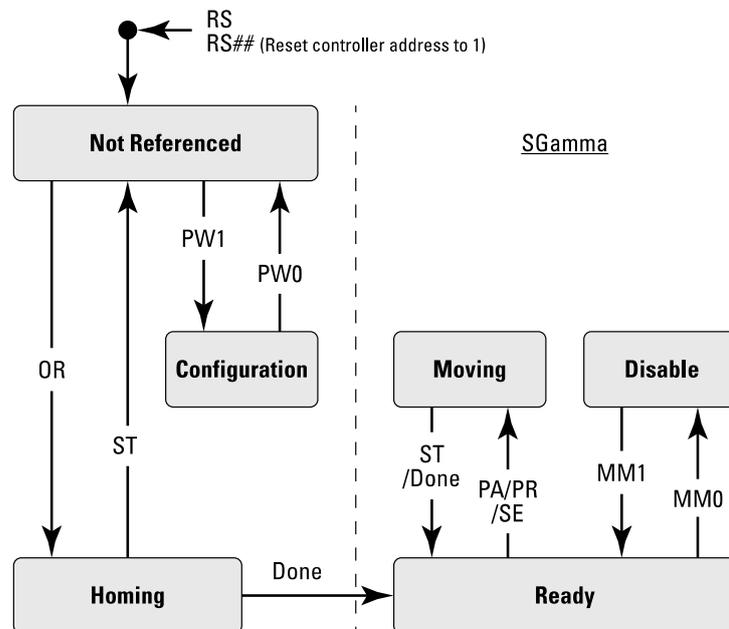

Actions in each state when End of Runs is encountered

NOT REFERENCED: No action.

CONFIGURATION: No action.

HOMING: Only check at end of HOMING and then change to NOT REFERENCED state.

MOVING: Abort motion and then changes to NOT REFERENCED state.

READY: Changes to NOT REFERENCED state.

DISABLE: Changes to NOT REFERENCED state.

After connecting the FC series stage to a power source, the stage is in the NOT REFERENCED state and must be initialized first. When initialization is successful, the controller goes to the NOT REFERENCED state. From the NOT REFERENCED state, the controller can go to the CONFIGURATION state with the PW1 command. In the CONFIGURATION state, the FC series stage allows changing stage and motor configuration parameters. The PW0 command saves all changes to the controller's memory and returns the controller back to the NOT REFERENCED state.

In the READY state, the motor is energized and ready to move. During a move execution (PA/PR), the controller is in the MOVING state and goes automatically back to the READY state when the move is completed. When errors are generated in the READY state, it changes the controller to the NOT REFERENCED state.

In the DISABLE state, the motor is not energized and move commands cannot be accepted. To go from the READY state to the DISABLE state and vice versa, use the MM command.

To go from the READY state or the DISABLE state back to the NOT REFERENCED state, for instance to make further parameter changes in the CONFIGURATION state, you must reset the controller with the RS command.



CAUTION

Changing stage and motor configuration parameters in the CONFIGURATION state is recommended only for knowledgeable person who understands the effect of the changes. See PW command for details.

4.2 Command Syntax

The FC series is a command-driven controller. The general format of a command is a two-letter ASCII word preceded and followed by parameters specific to the command:

Command format

nn	AA	xx
-----------	-----------	-----------

nn — Controller address, or
nothing if the issued command addresses all controllers.

AA — Command name.

xx — Parameter value, or
“?” to query the current value, or
nothing if the command takes no parameter.

Both upper and lower case characters are accepted. Depending on the command, it can have an optional or required prefix (**nn**) for the controller address and/or a suffix (**xx**) value, a “?” or no suffix at all.

Blank spaces

Blanks are allowed and ignored in any position, including inside a numerical value, unless enclosed within quotes. The following two commands are equivalent, but the first example might be confusing:

2P A1.43 6

2PA1.436

Decimal separator

A dot (“.”) must be used as decimal separator for all numerical values.

Command terminator

Commands are executed as either of the command terminator C_R or L_F (carriage-return, ASCII 13 or line-feed, ASCII 10) is received. The controller will analyze the received string. If the command is valid and its parameters are in the specified range, it will be executed. Otherwise it will memorize an error.

After the command parameters are identified, all remaining characters in the input string until the first command terminator, if any, will be ignored. Commands from the PC to the FC series may still be concatenated in a single string, but each command must be separated from the next one by a carriage-return or a line-feed.

In case any error occurs, the reported error will be recorded and can be checked using the TE command. Please refer to the command set in section 4.4 for details.

4.3 Command Execution Time

The FC series controller interprets commands continuously as they are received. The typical execution time for a "tell position command" (nTP?) is about 10 ms. Here, command execution time means the between sending a command and receiving an answer.

It is important to note that a move command, that may lasts for several seconds, will not suspend the controller from further command execution. So for an efficient process flow with many move commands it is recommended to use the PT command (get time for a relative move), and to query the controller status (TS command) or the current position (TP command) before any further motion command is sent.

4.4 Command Set

This section describes the supported two-letter ASCII commands used to configure and operate the FC series. The general command format is:

Command format

nn	AA	xx
-----------	-----------	-----------

nn — Optional or required controller address.

AA — Command name.

xx — Optional or required value or “?” to query current value.

Most commands can be used to set a value (in that case the command name is followed by the value, represented here as “**xx**”) or to query the current value (in that case the command name is followed by a “?”). When querying a value, the controller responds with the command it received followed by the queried value.

Examples:

1VA10 sets the velocity of the controller #1 to 10 units/second (and sends nothing back).

1VA? sends back the reply "1VA10", which means: "controller #1's velocity is 10 units/second".

Not every command can be executed in all states of the FC series and some commands have different meaning in different states. It is therefore important to understand the state diagram of the controller, see section 4.1.

	Not Ref.	Config.	Disable	Ready	Motion	Description
AC	–	○	□	□	–	Set/Get acceleration
BA	–	○	–	–	–	Set/Get backlash compensation
BH	–	○	–	–	–	Set/Get hysteresis compensation
FR	–	○	–	–	–	Set/Get stepper motor configuration
HT	–	○	–	–	–	Set/Get HOME search type
ID	–	○	□	□	–	Set/Get stage identifier
JR	–	○	□	□	–	Set/Get jerk time
MM	–	–	●	●	–	Enter/Leave DISABLE state
OH	–	○	–	–	–	Set/Get HOME search velocity
OR	●	–	–	–	–	Execute HOME search
OT	–	○	–	–	–	Set/Get HOME search time-out
PA	–	–	–	●	–	Move absolute
PR	–	–	–	●	–	Move relative
PT	–	–	●	●	●	Get estimated duration of a relative move
PW	●	●	–	–	–	Enter/Leave CONFIGURATION state
RS	●	●	●	●	●	Reset controller
RS##	●	●	●	●	●	Reset controller's address to 1
SA	–	○	–	–	–	Set/Get controller's RS-485 address
SE	–	–	–	●	–	Configure/Execute simultaneous started move
SL	–	○	□	□	–	Set/Get negative software limit
SR	–	○	□	□	–	Set/Get positive software limit
ST	–	–	–	–	●	Stop motion
TB	●	●	●	●	●	Get command error string
TE	●	●	●	●	●	Get last command error
TH	●	●	●	●	●	Get set-point position
TP	●	●	●	●	●	Get current position
TS	●	●	●	●	●	Get positioner error and controller state
VA	–	○	□	□	–	Set/Get velocity
VE	●	●	●	●	●	Get controller revision information
ZT	●	●	●	●	●	Get all axis parameters

Not Ref. Corresponds to the NOT REFERENCED state (for details see state diagram, section 4.1).

Config. Corresponds to the CONFIGURATION state.

Disable Corresponds to the DISABLE state.

Ready Corresponds to the READY state.

Motion Corresponds to the HOMING and MOVING states.

○ Changes configuration parameters. Those changes will be stored in the controller's memory with the PW1 command and remain available after switching off the controller.

□ Changes working parameters only. Those changes will get lost when switching off the controller.

● Accepted command.

– Command is forbidden in this state (will memorize an error).

Grey line Command passed without preceding controller number applies to all controllers (e.g. MM0 disables all controllers).

AC — Set/Get acceleration

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	□	□	–
Syntax	xxACnn or xxAC?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Acceleration value.			
Range	xx	—	1 to 31			
	nn	—	> 10⁻⁶ and < 10¹²			
Units	xx	—	None			
	nn	—	Preset units/s ²			
Defaults	xx Missing:		Error B.			
	Out of range:		Error B.			
	nn Missing:		Error C.			
	Out of range:		Error C.			
Description	<p>In CONFIGURATION state, this command sets the maximum acceleration value which can then be saved in the controller's non-volatile memory using the PW command. This is the maximum acceleration that can be applied to the mechanical system. It is also the default acceleration that will be used for all moves unless a lower value is set in DISABLE or READY state.</p> <p>In DISABLE or READY state, this command sets the acceleration used for all subsequent moves. Its value can be up to the programmed value in CONFIGURATION state. This value is not saved in the controller's memory and will be lost after reboot.</p>					
Returns	If the sign “?” is used instead of nn , this command returns the current value for the state in which the controller is (either CONFIGURATION or DISABLE/READY).					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	JR	—	Set/Get jerk time.			
	VA	—	Set/Get velocity.			
Example	1AC500		<i>Set controller #1 acceleration to 500 units/s².</i>			
	1AC?		<i>Controller returns 1AC500.</i>			

BA — Set/Get backlash compensation

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	–	–	–
Syntax	xxBAnn or xxBA?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Backlash value.			
Range	xx	—	1 to 31			
	nn	—	≥ 0 and $< 10^{12}$			
Units	xx	—	None			
	nn	—	Preset units			
Defaults	xx Missing:		Error B.			
	Out of range:		Error B.			
	nn Missing:		Error C.			
	Out of range:		Error C.			
Description	<p>The BA command sets the backlash compensation value. This is the value that the controller moves the motor in addition to the commanded distance with any move that reverses the direction of motion without changing the current position value (TP command).</p> <p>The BA command helps compensating for repeatable mechanical defects that appear when reversing the direction of motion, for instance mechanical wear. The value 0 disables this function. This feature can be only used when the hysteresis compensation (BH) is disabled.</p> <p>When a value different from 0 is set, the travel range of the stage is decreased by the same amount.</p>					
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	J	—	Execution not allowed in DISABLE state.			
	K	—	Execution not allowed in READY state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	BH	—	Set hysteresis compensation.			
Example	1BA0.005		<i>Set controller #1 backlash compensation to 0.005 units.</i>			

BH — Set/Get hysteresis compensation

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		—	○	—	—	—
Syntax	xxBHnn or xxBH?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Hysteresis value.			
Range	xx	—	1 to 31			
	nn	—	≥ 0 and $< 10^{12}$			
Units	xx	—	None			
	nn	—	Preset units			
Defaults	xx Missing:	Error B.				
	Out of range:	Error B.				
	Floating point:	Error A.				
	nn Missing:	Error C.				
	Out of range:	Error C.				
Description	<p>The BH command sets the hysteresis compensation value. When set to a value different than zero, the controller will issue for each move in the positive direction a move of the commanded distance plus the hysteresis compensation value, and then a second move of the hysteresis compensation value in the negative direction. This motion ensures that a final position gets always approached from the same direction and distance and helps compensating for non-repeatable mechanical defects like hysteresis or mechanical stiffness variations.</p> <p>The value 0 disables this function. The BH command can not be used when the backlash compensation is enabled (BA command).</p> <p>When a value different from 0 is set, the travel range of the stage is decreased by the same amount in the positive direction.</p>					
NOTE						
<p>The homing set on the positive end of run and hysteresis compensation are not compatible. Any attempt to use both features together will make the stage fail.</p>						
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	J	—	Execution not allowed in DISABLE state.			
	K	—	Execution not allowed in READY state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	BA	—	Set backlash compensation.			
Example	1BH0.015		<i>Set controller #1 backlash compensation to 0.015 units.</i>			

FR — Set/Get stepper motor configuration

Usage	Not Ref.	Config.	Disable	Ready	Motion
	–	○	–	–	–
Syntax	xxFRSnn, xxFRM? or xxFRS?				
Parameters					
Description	xx [int]	—	Controller address.		
	Mmm [int]	—	Amount of micro-steps per full step.		
	Snn [float]	—	Full step displacement length in 1/1000 of unit.		
Range	xx	—	1 to 31		
	mm	—	> 0 and ≤ 2000		
	nn	—	> 10⁻⁶ and < 10¹²		
Units	xx	—	None.		
	Mmm	—	None.		
	Snn	—	1/1000 of unit.		
Defaults	xx Missing:	Error B.			
	Out of range:	Error B.			
	mm Missing:	Error C.			
	Out of range:	Error C.			
	nn Missing:	Error C.			
	Out of range:	Error C.			
Description	FRM: For compatibility. No effect. Always 128 µsteps.				
	FRS: This command sets the displacement length per full step in 1/1000 of unit.				
Returns	If the sign “?” is used instead of mm or nn , this command returns the current programmed value.				
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
	C	—	Parameter missing or out of range.		
	D	—	Execution not allowed.		
	H	—	Execution not allowed in NOT REFERENCED state.		
	J	—	Execution not allowed in DISABLE state.		
	K	—	Execution not allowed in READY state.		
	L	—	Execution not allowed in HOMING state.		
	M	—	Execution not allowed in MOVING state.		
Example	1FRS10		<i>Set controller #1 full step value to 10 milli-units.</i>		

HT — Set/Get HOME search type

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	–	–	–
Syntax	xxHTnn or xxHT?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [int]	—	Home search type identifier.			
Range	xx	—	1 to 31			
	nn	—	1 use current position as HOME. 2 use MZ switch (mechanical zero) to detect HOME position. 4 use EoR- switch (negative end of range) to detect HOME position.			
Units	xx	—	None.			
	nn	—	None.			
Defaults	xx Missing:		Error B.			
	Out of range:		Error B.			
	nn Missing:		Error C.			
	Out of range:		Error C.			
Description	This command sets the type of HOME search used with the OR command.					
NOTE						
The homing set on the positive end of run and hysteresis compensation are not compatible. Any attempt to use both features together will make the stage fail.						
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	J	—	Execution not allowed in DISABLE state.			
	K	—	Execution not allowed in READY state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	OH	—	Set/Get HOME search velocity.			
	OR	—	Execute HOME search.			
	OT	—	Set HOME search time-out.			
Example	1HT1		<i>Set controller #1 HOME sequence to use current position.</i>			

ID — Set/Get stage identifier

Usage	Not Ref.	Config.	Disable	Ready	Motion
	–	○	□	□	–
Syntax	xxIDnn or xxID?				
Parameters					
Description	xx [int]	—	Controller address.		
	nn [string]	—	Stage identifier string.		
Range	xx	—	1 to 31		
	nn	—	1 to 31 ASCII characters.		
Units	xx	—	None		
	nn	—	None		
Defaults	xx Missing:	Error B.			
	Out of range:	Error B.			
	nn Missing:	Error C.			
	Out of range:	Error C.			
Description	<p>The ID command sets the stage identifier in the form of a character string. Any printable character can be used ; spaces are admissible only if the string is enclosed in quotes, "like this". If not, spaces and tabs are ignored.</p> <p>In CONFIGURATION state, this command sets a new value for the stage identifier which can then be saved in the device's non-volatile memory with the PW command. It is also the default value that will be used unless a different value is set in DISABLE or READY state.</p> <p>In DISABLE or READY state, this command allows setting a new working parameter for the stage identifier. This value is not saved in the controller's memory and will be lost after reboot.</p>				
Returns	<p>If the sign "?" is used instead of nn, this command returns the current identifier string for the state in which the controller is (either CONFIGURATION or DISABLE/READY).</p>				
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
	C	—	Parameter missing or out of range.		
	D	—	Execution not allowed.		
	H	—	Execution not allowed in NOT REFERENCED state.		
	L	—	Execution not allowed in HOMING state.		
	M	—	Execution not allowed in MOVING state.		
Example	<pre>IID? Get stage identifier for controller #1. IID URS100CC Set controller #1's stage identifier to: URS100CC.</pre>				

JR — Set/Get jerk time

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	□	□	–
Syntax	xxJRnn or xxJR?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Jerk time value.			
Range	xx	—	1 to 31			
	nn	—	> 0.001 and < 10¹²			
Units	xx	—	None.			
	nn	—	Seconds.			
Defaults	xx Missing:		Error B.			
	Out of range:		Error B.			
	nn Missing:		Error C.			
	Out of range:		Error C.			
Description	<p>Jerk is the derivative of acceleration. The jerk time defines the time to reach the needed acceleration. A longer jerk time reduces stress to the mechanics and smoothes motion.</p> <p>In CONFIGURATION state, this command sets the value for the maximum jerk time which can then be saved in the controller's non-volatile memory using the PW command. It is also the default value that will be used unless a different value is set in DISABLE or READY state.</p> <p>In DISABLE or READY state, this command allows setting a new working parameter for the maximum jerk time. This value is not saved in the controller's memory and will be lost after reboot.</p>					
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value for the state in which the controller is (either CONFIGURATION or DISABLE/READY).					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution impossible (axis in movement).			
	H	—	Execution not allowed in NOT REFERENCED state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	AC	—	Set/Get acceleration.			
	VA	—	Set/Get velocity.			
Example	1JR0.05		<i>Set controller #1 jerk time to 0.05 seconds.</i>			

MM — Enter/Leave DISABLE state

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		—	—	●	●	—
Syntax	xxMMnn or xxMM?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [int]	—	Whether to enter (1) or leave (0) the DISABLE state.			
Range	xx	—	0 to 31			
	nn	—	0 changes state from READY to DISABLE. 1 changes state from DISABLE to READY.			
Units	xx	—	None.			
	nn	—	None.			
Defaults	xx	Missing:	Change to 0 (will forward this command to all controllers).			
		Out of range:	Error B.			
	nn	Missing:	Error C.			
		Out of range:	Error C.			
Description	When the MM command is sent without preceding controller number or the controller number is 0, the MM command is executed on all controllers.					
	MM0 changes the controller's state from READY to DISABLE. In DISABLE state the control loop is open and the motor is not powered.					
	MM1 changes the controller's state from DISABLE to READY. The controller's set point position is set equal to its current position and the control loop gets closed (depending on the closed-loop state). Any residual error is cleared and the motor is powered.					
Returns	If the sign “?” is used instead of nn , this command returns the current state. Refer to the TS command section for the list of controller states.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	I	—	Execution not allowed in CONFIGURATION state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	PW	—	Enter/leave CONFIGURATION state.			
Example	MM0		<i>All controllers go to DISABLE state.</i>			

OH — Set/Get HOME search velocity

Usage	Not Ref.	Config.	Disable	Ready	Motion
	–	○	–	–	–
Syntax	xxOHnn or xxOH?				
Parameters					
Description	xx [int]	—	Controller address.		
	nn [float]	—	HOME search velocity.		
Range	xx	—	1 to 31		
	nn	—	> 10⁻⁶ and < 10¹²		
Units	xx	—	None.		
	nn	—	Preset units/s.		
Defaults	xx Missing:	Error B.			
	Out of range:	Error B.			
	nn Missing:	Error C.			
	Out of range:	Error C.			
Description	This command sets the maximum velocity used by the controller for the HOME search.				
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value.				
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
	C	—	Parameter missing or out of range.		
	D	—	Execution not allowed.		
	H	—	Execution not allowed in NOT REFERENCED state.		
	J	—	Execution not allowed in DISABLE state.		
	K	—	Execution not allowed in READY state.		
	L	—	Execution not allowed in HOMING state.		
	M	—	Execution not allowed in MOVING state.		
Rel. Commands	HT	—	Set/Get HOME search type.		
	OR	—	Execute HOME search.		
	OT	—	Set HOME search time-out.		
Example	1OH50		<i>Set controller #1 HOME search velocity to 50 units/s.</i>		

OR — Execute HOME search

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		●	—	—	—	—
Syntax	xxOR					
Parameters						
Description	xx [int]	—	Controller address.			
Range	xx	—	1 to 31			
Units	xx	—	None.			
Defaults	xx Missing:		Error B.			
	Out of range:		Error B.			
	nn Missing:		Error C.			
	Out of range:		Error C.			
Description	This command starts the execution of the HOME search according to the algorithm defined by the HT command.					
	When in NOT REFERENCED state, for instance after system start, any positioner must first be homed with the OR command before further motion commands can be executed.					
	The OR command is accepted only in NOT REFERENCED state and only when no hardware error is present (except end-of-runs). Refer to the TS command to get more information on the possible hardware errors.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	E	—	home sequence already started.			
	I	—	Execution not allowed in CONFIGURATION state.			
	J	—	Execution not allowed in DISABLE state.			
	K	—	Execution not allowed in READY state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	HT	—	Set HOME search type.			
	OH	—	Set HOME search velocity.			
	OT	—	Set HOME search time-out.			
Example	1OR		<i>Execute HOME search with controller #1.</i>			

OT — Set/Get HOME search time-out

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	–	–	–
Syntax	xxOTnn or xxOT?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	HOME time-out.			
Range	xx	—	1 to 31			
	nn	—	> 1 and < 1000			
Units	xx	—	None.			
	nn	—	Seconds			
Defaults	xx Missing:		Error B.			
	Out of range:		Error B.			
	nn Missing:		Error C.			
	Out of range:		Error C.			
Description	This command sets the time-out value for the HOME search. When the HOME search does not finish successfully before this delay elapses, the HOME search is aborted and an error is recorded.					
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	J	—	Execution not allowed in DISABLE state.			
	K	—	Execution not allowed in READY state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	HT	—	Set HOME search type.			
	OH	—	Set HOME search velocity.			
	OR	—	Execute HOME search.			
Example	1OT2.2		<i>Set controller #1 HOME time-out to 2.2 seconds.</i>			

PA — Move absolute

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		—	—	—	●	—
Syntax	xxPAnn or xxPA?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	New absolute position.			
Range	xx	—	1 to 31			
	nn	—	\geq SL and \leq SR			
Units	xx	—	None.			
	nn	—	Preset units.			
Defaults	xx Missing:	Error B.				
	Out of range:	Error B.				
	nn Missing:	Error C.				
	Out of range:	Error C.				
Description	The PA command initiates an absolute move. When received, the positioner will move, with the predefined acceleration and velocity, to the new absolute position specified by nn .					
	The PA command is only accepted in READY state, AND when the new absolute position is higher or equal to the negative software limit (SL), AND lower or equal to the positive software limit (SR).					
	The controller always rounds the new target position to the closest micro-step position.					
Returns	If the sign “?” is used instead of nn , this command returns the target absolute position value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	G	—	Target position out of limits.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	I	—	Execution not allowed in CONFIGURATION state.			
	J	—	Execution not allowed in DISABLE state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	PR	—	Move relative.			
	PT	—	Get motion time for a relative move.			
	TH	—	Get set-point position.			
	TP	—	Get current position.			
Example	1PA2.2		<i>Move positioner on controller #1 to absolute position 2.2 units.</i>			

PR — Move relative

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		—	—	—	●	—
Syntax	xxPRnn or xxPR?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Displacement.			
Range	xx	—	1 to 31			
	nn	—	$\geq (\text{SL} - \text{TP})$ and $\leq (\text{SR} - \text{TP})$			
Units	xx	—	None.			
	nn	—	Preset units.			
Defaults	xx Missing:	Error B.				
	Out of range:	Error B.				
	nn Missing:	Error C.				
	Out of range:	Error C.				
Description	The PR command initiates a relative move. When received, the positioner will move, with the predefined acceleration and velocity, to a new absolute position nn units away from the current absolute position.					
	The PR command gets only accepted in READY state, AND when the distance of the positioner to the software limit in the same direction is longer than the commanded displacement.					
	The controller always rounds the new target position to the closest micro-step position.					
Returns	If the sign “?” is used instead of nn , this command returns the target absolute position value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	G	—	Displacement out of limits.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	I	—	Execution not allowed in CONFIGURATION state.			
	J	—	Execution not allowed in DISABLE state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	PA	—	Move absolute.			
	PT	—	Get motion time for a relative move.			
	TH	—	Get set-point position.			
	TP	—	Get current position.			
Example	1PR2.2		<i>Move positioner on controller #1 to a new position 2.2 units away from the current position.</i>			

PT — Get motion time for a relative move

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		—	—	●	●	●
Syntax	xxPTnn					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Displacement.			
Range	xx	—	1 to 31			
	nn	—	> 10⁻⁶ and < 10¹²			
Units	xx	—	None.			
	nn	—	Preset units.			
Defaults	xx Missing:	Error B.				
	Out of range:	Error B.				
	nn Missing:	Error C.				
	Out of range:	Error C.				
Description	The PT commands helps evaluating move times for an efficient program flow. When receiving the PT command, the controller computes and returns the time, in seconds, necessary to execute a relative move of displacement nn with the current working parameters (velocity, acceleration, etc.). The controller does not execute any displacement.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	I	—	Execution not allowed in CONFIGURATION state.			
Rel. Commands	PA	—	Move absolute.			
	PR	—	Move relative.			
	TH	—	Get set-point position.			
	TP	—	Get current position.			
Example	1PT2.2		<i>Get time to move positioner on controller #1 by 2.2 units.</i>			
	1PT0.25		<i>Controller returns: 0.25 seconds.</i>			

PW — Enter/Leave CONFIGURATION state

Usage	Not Ref.	Config.	Disable	Ready	Motion
	●	●	—	—	—
Syntax	xxPWnn or xxPW?				
Parameters					
Description	xx [int]	—	Controller address.		
	nn [int]	—	Whether to enter (1) or leave (0) the CONFIGURATION state.		
Range	xx	—	1 to 31		
	nn	—	1: Go from NOT REFERENCED state to CONFIGURATION state. 0: Go from CONFIGURATION state to NOT REFERENCED state.		
Units	xx	—	None.		
	nn	—	None.		
Defaults	xx	Missing:	Error B.		
		Out of range:	Error B.		
	nn	Missing:	Error C.		
		Out of range:	Error C.		
Description	<p>PW1 changes the controller's state from NOT REFERENCED to CONFIGURATION. In CONFIGURATION state, all parameter settings are saved in the controller's memory upon exiting this state and remain available after switching off the controller. In addition, some settings are only possible in CONFIGURATION state (e.g. set drive voltage, set Backlash compensation, etc.).</p> <p>PW0 checks all stage parameters, and if they are acceptable, saves them in the flash memory of the controller. After that, it changes the controller's state from CONFIGURATION to NOT REFERENCED.</p> <p>The execution of a PW0 command may take up to 5 seconds. During that time the controller will not respond to any other command.</p>				
Returns	If the sign "?" is used instead of nn, this command returns whether we are or not in the CONFIGURATION state.				
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
	C	—	Parameter missing or out of range.		
	D	—	Execution not allowed.		
	J	—	Execution not allowed in DISABLE state.		
	K	—	Execution not allowed in READY state.		
	L	—	Execution not allowed in HOMING state.		
	M	—	Execution not allowed in MOVING state.		
Rel. Commands	MM	—	Enter/Leave DISABLE state.		
Example	1PW1		<i>Changes controller #1 to CONFIGURATION state.</i>		

NOTE

The PW command is limited to 100 writes. Unit failure due to excessive use of the PW command is not covered by the warranty.

The PW command is used to change the default configuration parameters that are stored in memory, and not working parameters that may be changed on the fly.

RS — Reset controller

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		●	●	●	●	●
Syntax	xxRS					
Parameters						
Description	xx [int]	—	Controller address.			
Range	xx	—	1 to 31			
Units	xx	—	None.			
Defaults	xx	Missing:	Error B.			
		Out of range:	Error B.			
Description	The RS command issues a hardware reset of the controller, equivalent to a power cycle. To go from DISABLE or READY state to CONFIGURATION state, it is also needed to reset the controller with the RS command, and then to change the controller's state with the PW1 command from NOT REFERENCED to CONFIGURATION.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	D	—	Execution not allowed.			
Example	1RS		<i>Reset controller #1.</i>			

RS## — Reset controller's address

Usage	Not Ref.	Config.	Disable	Ready	Motion
	○	○	○	○	○
Syntax	xxRS## or RS##				
Parameters					
Description	xx [int]	—	Controller address.		
Range	xx	—	0 to 31		
Units	xx	—	None.		
Defaults	xx	Missing:	Change to 0 (will forward this command to all controllers).		
		Out of range:	Error B.		
Description	The RS## command resets the controller's address to 1. This address needs to be different for each FC series when connected on a RS-485 communication network.				
Returns					
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
	D	—	Execution not allowed.		
Rel. Commands	SA	—	Set/Get controller's RS-485 address.		
Example	RS##		<i>Reset all controllers' addresses to 1.</i>		

SA — Set/Get controller's RS422 address

Usage	Not Ref.	Config.	Disable	Ready	Motion
	–	○	–	–	–
Syntax	xxSAnn or xxSA?				
Parameters					
Description	xx [int]	—	Current controller address.		
	nn [int]	—	New controller address (a.k.a. RS422 address).		
Range	xx	—	1		
	nn	—	1 to 31		
Units	xx	—	None.		
	nn	—	None.		
Defaults	xx Missing:	Error B.			
	Out of range:	Error B.			
	nn Missing:	Error C.			
	Out of range:	Error C.			
Description	The SA command sets the controller's RS422 address, also known as the axis number. This address needs to be different for each FC series when connected on a RS422 communication network.				
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value.				
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
	C	—	Parameter missing or out of range.		
	D	—	Execution not allowed.		
	H	—	Execution not allowed in NOT REFERENCED state.		
	J	—	Execution not allowed in DISABLE state.		
	K	—	Execution not allowed in READY state.		
	L	—	Execution not allowed in HOMING state.		
	M	—	Execution not allowed in MOVING state.		
Rel. Commands	RS##	—	Reset controller's address.		
Example	1SA3		<i>Set controller's RS422 address to 3.</i>		

SE — Configure/Execute simultaneous started move

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		—	—	—	●	—
Syntax	xxSEnn, xxSE? or SE					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	New target position.			
Range	xx	—	0 to 31			
	nn	—	≥ SL and ≤ SR			
Units	xx	—	None.			
	nn	—	Preset units.			
Defaults	xx	Missing:	Change to 0 (will forward this command to all controllers).			
		Out of range:	Error B.			
	nn	Missing:	Error C.			
		Out of range:	Error C.			
Description	<p>The SE command allows starting a move on different controllers at the same time.</p> <p>The command xxSEnn sets a new target position for the controller nn. But contrarily to the PA/PR commands, the move is not executed immediately, but only after receipt of an SE command without preceding controller number nor position value. When receiving the 2nd SE command, all controllers start moving to their target position.</p> <p>The xxSEnn command is only accepted in READY state, AND when the new target position is higher or equal to the negative software limit (SL), AND lower or equal to the positive software limit (SR). The controller always rounds the new target position to the closest micro-step position.</p> <p>The SE command should not be confused with a synchronized move. With a synchronized move, all positioners start moving simultaneously AND complete their moves at the same time. The SE command starts a move on all controllers at the same time, but each positioner moves with its individually defined velocity and acceleration. Consequently, the different positioners do NOT complete their moves at the same time.</p>					
Returns	If the sign “?” is used instead of nn, this command returns the target position value set by the SE command, which may not be the target position set by the PA/PR commands.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	I	—	Execution not allowed in CONFIGURATION state.			
	J	—	Execution not allowed in DISABLE state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	PA / PR	—	Move absolute / relative.			
	TH	—	Get set-point position.			
	TP	—	Get current position.			
Example	1SE2.2		<i>Prepare controller #1 to move to absolute position 2.2 units.</i>			
	2SE3.3		<i>Prepare controller #2 to move to absolute position 3.3 units.</i>			
	SE		<i>All controllers start their programmed move, if any.</i>			

SL — Set/Get negative software limit

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	□	□	–
Syntax	xxSLnn or xxSL?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Negative software limit.			
Range	xx	—	1 to 31			
	nn	—	$> -10^{12}$ and ≤ 0			
Units	xx	—	None.			
	nn	—	Preset units.			
Defaults	xx	Missing:	Error B.			
		Out of range:	Error B.			
	nn	Missing:	Error C.			
		Out of range:	Error C.			
Description	<p>In CONFIGURATION state, this command sets the negative software limit which can then be saved in the controller's non-volatile memory using the PW command. It is also the default value that will be used unless a different value is set in DISABLE or READY state.</p> <p>In DISABLE or READY state, this command allows setting a new working parameter for the negative software limit. It must be lower than or equal to the set-point position. This value is not saved in the controller's memory and will be lost after reboot.</p> <p>The software limits are used to restrict the travel range of a positioner. For FC series, it is not possible to disable software limits. For an almost infinite motion, for instance with a FCR100 stage, set the lowest possible value, which is -99999930400.</p>					
Returns	If the sign “?” is used instead of nn, this command returns the current programmed value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	SR	—	Set positive software limit.			
Example	1SL-100		<i>Set controller #1 negative software limit to -100 units.</i>			

NOTE

For FCR100, the maximum allowed values for software limits are $\pm 2 \times 10^{14}$.

SR — Set/Get positive software limit

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	□	□	–
Syntax	xxSRnn or xxSR?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Positive software limit.			
Range	xx	—	1 to 31			
	nn	—	≥ 0 and $< 10^{12}$			
Units	xx	—	None.			
	nn	—	Preset units.			
Defaults	xx	Missing:	Error B.			
		Out of range:	Error B.			
	nn	Missing:	Error C.			
		Out of range:	Error C.			
Description	<p>In CONFIGURATION state, this command sets the positive software limit which can than be saved in the controller's non-volatile memory using the PW command. It is also the default value that will be used unless a different value is set in DISABLE or READY state.</p> <p>In DISABLE or READY state, this command allows setting a new working parameter for the positive software limit. It must be greater or equal to the set-point position. This value is not saved in the controller's memory and will be lost after reboot.</p> <p>The software limits are used to restrict the travel range of a positioner. For FC series, it is not possible to disable software limits. For an almost infinite motion, for instance with a FCR100 stage, set the highest possible value, which is 999999930400.</p>					
Returns	If the sign “?” is used instead of nn , this command returns the current programmed value.					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	SL	—	Set negative software limit.			
Example	1SR100		<i>Set controller #1 positive software positive to 100 units.</i>			

NOTE

For FCR100, the maximum allowed values for software limits are $\pm 2 \times 10^{14}$.

ST — Stop motion

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		—	—	—	—	●
Syntax	[xx]ST					
Parameters						
Description	xx [int]	—	Controller address.			
Range	xx	—	0 to 31			
Units	xx	—	None.			
Defaults	xx	Missing:	Change to 0 (will forward this command to all controllers).			
		Out of range:	Error B.			
Description	<p>The ST command is a safety feature. It stops a move in progress by decelerating the positioner immediately with the acceleration defined by the AC command until it stops.</p> <p>The xxST command with a preceding controller address stops a move in progress on controller xx. The ST command without any preceding controller address stops the moves on ALL controllers.</p>					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	I	—	Execution not allowed in CONFIGURATION state.			
	J	—	Execution not allowed in DISABLE state.			
	K	—	Execution not allowed in READY state.			
Example	ST		<i>Stop moves on all controllers.</i>			

TB — Get command error string

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		●	●	●	●	●
Syntax	xxTBnn					
Parameters						
Description	xx [int]	—	Controller address.			
Range	xx	—	1 to 31			
	nn [char]	—	Error code (refer to TE command).			
Units	xx	—	None.			
Defaults	xx	Missing:	Error B.			
		Out of range:	Error B.			
	nn	Missing:	Returns explanation of current error as a literal string.			
		Out of range:	Error C.			
Description	The TB command returns a literal character string that explains the meaning of the error code nn (see TE command for complete list).					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
Rel. Commands	TE	—	Get last command error.			

Example 1TB@ | *Get explanation to error code @.*
 1TB@ No error | *Controller returns: @ = means no error.*

TE — Get last command error

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		●	●	●	●	●
Syntax	xxTE					
Parameters						
Description	xx [int]	—	Controller address.			
Range	xx	—	1 to 31			
Units	xx	—	None.			
Defaults	xx	Missing:	Error B.			
		Out of range:	Error B.			
Description	<p>The TE command returns the currently recorded error. When a command is not executable or fails, an error is recorded. This error can be read with the TE command. After the execution of a TE command, the error buffer is erased and another TE command will return @, which means "No error". When a new command error is generated before the previous command error is read, the new command's error will overwrite the currently memorized error.</p> <p>For a safe program flow it is recommended to always query the command error after each command execution.</p>					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	D	—	Execution not allowed.			
Rel. Commands	TB	—	Get command error string.			
Example	1TE		<i>Get last error memorized on controller #1.</i>			
			<i>Controller returns: 1TE@, means no error.</i>			

List of errors and corresponding strings (see TB command):

@	—	No error.
A	—	Unknown message code or floating point controller address.
B	—	Controller address not correct.
C	—	Parameter missing or out of range.
D	—	Command not allowed.
E	—	Home sequence already started.
G	—	Displacement out of limits.
H	—	Command not allowed in NOT REFERENCED state.
I	—	Command not allowed in CONFIGURATION state.
J	—	Command not allowed in DISABLE state.
K	—	Command not allowed in READY state.
L	—	Command not allowed in HOMING state.
M	—	Command not allowed in MOVING state.
N	—	Current position out of software limit.
S	—	Communication Time Out.
U	—	Error during EEPROM access.
V	—	Error during command execution.

TH — Get set-point position

Usage	Not Ref.	Config.	Disable	Ready	Motion
	●	●	●	●	●
Syntax	xxTH				
Parameters					
Description	xx [int]	—	Controller address.		
Range	xx	—	1 to 31		
Units	xx	—	None.		
Defaults	xx	Missing:	Error B.		
		Out of range:	Error B.		
Description	The TH command returns the value of the set-point or theoretical position. This is the position where the positioner should be. In MOVING state, the set-point position changes according to the calculation of the motion profiler. In READY state, the set-point position is equal to the target position.				
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
	D	—	Execution not allowed.		
Rel. Commands	TP	—	Get current position.		
Example	1TH		<i>Get set-point position of controller #1.</i>		
	1TH0		<i>Controller returns: set-point position = 0 units.</i>		

TP — Get current position

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		●	●	●	●	●
Syntax	xxTP					
Parameters						
Description	xx [int]	—	Controller address.			
Range	xx	—	1 to 31			
Units	xx	—	None.			
Defaults	xx	Missing:	Error B.			
		Out of range:	Error B.			
Description	<p>The TP command returns the value of the current position. This is the position where the positioner actually is. In MOVING state, this value always changes. In READY state, this value should be equal or very close to the set-point and target position.</p> <p>Together with the TS command, the TP command helps evaluating whether a motion has completed.</p>					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	D	—	Execution not allowed			
Rel. Commands	TH	—	Get set-point position.			
	TS	—	Get positioner error and controller state.			
Example	1TP		<i>Get current position of controller #1.</i>			
	1TP0		<i>Controller returns: actual position = 0 units.</i>			

TS — Get positioner error and controller state

Usage	Not Ref.	Config.	Disable	Ready	Motion
	●	●	●	●	●
Syntax	xxTS				
Parameters					
Description	xx [int]	—	Controller address.		
Range	xx	—	1 to 31		
Units	xx	—	None.		
	nn	—	None.		
Defaults	xx	Missing:	Error B.		
		Out of range:	Error B.		

Description The TS command returns the positioner error and the current controller state.

Returns The TS command returns six characters (1TSabcdef). The first 4 characters (abcd) represent the positioner error as hexadecimal number. The last two characters (ef) represent the controller state as hexadecimal number.

Error code (abcd): Convert each hexadecimal to a binary:

F	E	D	C	B	A	9	8	7	6	5	4	3	2	1	0
1111	1110	1101	1100	1011	1010	1001	1000	0111	0110	0101	0100	0011	0010	0001	0000

E

each bit represents one possible error (exception made of bit C1):

A				B				C				D			
1	1	1	1	1	1	1	1	1	1	1	1	1	1	1	1
Not used	Not used	Not used	Not used	Driver overheating	Driver fault	Not used	Not used	No parameters in memory	Homing time out	Not used	Newport reserved (MZ status)	RMS current limit	Not used	Positive end of run	Negative end of run

NOTE

Bit C1 (MZ status) is not an error. It is reserved for Newport technicians to diagnose the mechanical zero sensor status during customer support and servicing.

Examples:

- Error map 0000 = No errors
- Error map 0002 = Positive end of run
- Error map 0048 = Homing time out, RMS current limit

Controller states (ef):

- **0A**: NOT REFERENCED from RESET.
- **0B**: NOT REFERENCED from HOMING.
- **0C**: NOT REFERENCED from CONFIGURATION.
- **0D**: NOT REFERENCED from DISABLE.
- **0E**: NOT REFERENCED from READY.
- **0F**: NOT REFERENCED from MOVING.
- **10**: NOT REFERENCED - NO PARAMETERS IN MEMORY.
- **14**: CONFIGURATION.
- **1E**: HOMING.
- **28**: MOVING.
- **32**: READY from HOMING.
- **33**: READY from MOVING.
- **34**: READY from DISABLE.
- **3C**: DISABLE from READY.
- **3D**: DISABLE from MOVING.

NOTE

The positioner error gets updated periodically, approx. every 1 ms.

The TS command reads the positioner error and clears it at the same time (same as what the command TE does with command errors). So when launching the TS command, it is important to process the TS feedback accordingly.

Errors	A	—	Unknown message code or floating point controller address.
	B	—	Controller address not correct.
Rel. Commands	TE	—	Get last command error.

Example	1TS		<i>Get error and state of controller #1.</i>
	ITS00000A		<i>Controller returns: no errors and state is NOT REFERENCED from reset.</i>

VA — Set/Get velocity

	Usage	Not Ref.	Config.	Disable	Ready	Motion
		–	○	□	□	–
Syntax	xxVAnn or xxVA?					
Parameters						
Description	xx [int]	—	Controller address.			
	nn [float]	—	Velocity value.			
Range	xx	—	1 to 31			
	nn	—	> 10⁻⁶ and < 10¹²			
Units	xx	—	None.			
	nn	—	Preset units/s.			
Defaults	xx Missing:		Error B.			
	Out of range:		Error B.			
	nn Missing:		Error C.			
	Out of range:		Error C.			
Description	In CONFIGURATION state, this command sets the maximum (i.e. cruise) velocity value which can then be saved in the controller's non-volatile memory using the PW command. This should be the maximum velocity that can be applied to the mechanical system. It is also the default velocity that will be used for all moves unless a lower value is set in DISABLE or READY state.					
	In DISABLE or READY state, this command sets the velocity used for all subsequent moves. Its value can be up to the programmed value set in the CONFIGURATION state. This value is not saved in the controller's memory and will be lost after reboot.					
Returns	If the sign “?” is used instead of nn , this command returns the current value for the state in which the controller is (either CONFIGURATION or DISABLE/READY).					
Errors	A	—	Unknown message code or floating point controller address.			
	B	—	Controller address not correct.			
	C	—	Parameter missing or out of range.			
	D	—	Execution not allowed.			
	H	—	Execution not allowed in NOT REFERENCED state.			
	L	—	Execution not allowed in HOMING state.			
	M	—	Execution not allowed in MOVING state.			
Rel. Commands	AC	—	Set/Get acceleration.			
	JR	—	Set/Get jerk time.			
Example	1VA50		<i>Set controller #1 maximum velocity to 50 units/s.</i>			

VE — Get controller revision information

Usage	Not Ref.	Config.	Disable	Ready	Motion
	●	●	●	●	●
Syntax	xxVE				
Parameters					
Description	xx [int]	—	Controller address.		
	nn [string]	—	Action.		
Range	xx	—	1 to 31		
Units	xx	—	None.		
Defaults	xx	Missing:	Error B.		
		Out of range:	Error B.		
Description	This command returns the controller's firmware revision information.				
Errors	A	—	Unknown message code or floating point controller address.		
	B	—	Controller address not correct.		
Example	1VE		<i>Get controller #1 revision information.</i>		
	<i>1VE FC family controller 2.0.0 Controller returns revision number</i>				

ZT — Get all configuration parameters

Usage	Not Ref.	Config.	Disable	Ready	Motion
	●	●	●	●	●
Syntax	xxZT				
Parameters					
Description	xx [int]	—	Controller address.		
Range	xx	—	1 to 31		
Units	xx	—	None.		
Defaults	xx Missing:	Error B.			
	Out of range:	Error B.			
Description	The ZT command returns the list of all current configuration parameters. The ZT command allows a quick review of all current stage parameter and simplifies the configuration of Newport stages, for instance by copying all the returned values into a configuration file which can be later fed back to the stage by simply pasting its contents in the terminal emulator when the stage needs to be quickly reconfigured.				
Errors	A	—	Unknown message code or floating point controller address		
	B	—	Controller address not correct		
Example	1ZT		<i>Get controller #1 configuration data.</i>		
	1PW1				
	1AC320.000000				
	1BA0.000000				
	...				
	1VA80.000000				
	1PW0				

5.0 Maintenance

5.1 Maintenance

The FC series stages require no particular maintenance and there are no user-serviceable parts or user adjustments to be made. However, as with other precision mechanical stages, care must be taken for handling, operation and storage.



CAUTION

FCR100 stages must be used or stocked in a clean environment to avoid dust, humidity, solvents or other substances.

RECOMMENDATION

It is recommended to return your stage to Newport's After Sales Service after every 18 months of use for lubrication and inspection.

If your FC stage is mounted on a workstation and cannot be easily removed, please contact Newport's After Sales Service for further instructions.

5.2 Repair



CAUTION

Never attempt to disassemble any part of the stage unless it is described in this manual.

Disassembly of components, unless instructed by Newport, may lead to failure or malfunction of the stage.

If you observe any malfunction in your stage, please contact us immediately to arrange a repair service.



CAUTION

Any attempt to disassemble or repair the stage without authorization will void your warranty.

5.3 Troubleshooting

A list of the most common problems and their corrective actions is provided in the table below. Use this as a reference but remember that problems may be related to other operator errors and may be resolved by other simple solutions.

Problem (Error Message)	Cause	Corrective Actions
Please verify the stage is powered, then relaunch the user interface	Power cord not plugged in.	Plug the power cord in the appropriate outlet. Verify the 3-pin connector is seated properly.
	Power supply cable was connected to outlet before the USB cable was connected to PC.	The power supply cable must remain fully disconnected from the outlet when the USB cable is connected to the PC. Unplug both the USB cable and the power supply cable. Then, make the USB connection first.
	Bad connection	Check the integrity of cables and make sure each connector is fully secured.
Access to com port is denied.	Com port of the PC may be already used by another program.	Make sure to terminate all other programs (ie. LabVIEW driver) before launching the GUI.
Applet launching unsuccessful. Access to the log file denied.	When you are not the Administrator of the PC, Window 7 disables writing in the Program Files folder.	Obtain the full Admin right in the PC. Or simply copy the "Motion Control" folder from "Program Files" to "My Documents" folder.



Visit Newport Online at:
www.newport.com

North America & Asia

Newport Corporation
1791 Deere Ave.
Irvine, CA 92606, USA

Sales

Tel.: (800) 222-6440
e-mail: sales@newport.com

Technical Support

Tel.: (800) 222-6440
e-mail: tech@newport.com

Service, RMAs & Returns

Tel.: (800) 222-6440
e-mail: service@newport.com

Europe

MICRO-CONTROLE Spectra-Physics S.A.S
9, rue du Bois Sauvage
91055 Évry CEDEX
France

Sales

Tel.: +33 (0)1.60.91.68.68
e-mail: france@newport.com

Technical Support

e-mail: tech_europe@newport.com

Service & Returns

Tel.: +33 (0)2.38.40.51.55

